




FREEDOMX™

THE FREEDOMX™ COUCH EXTENSION MODULE

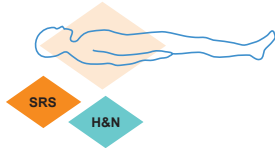

For an unobstructed 360° treatment of a patient's head and neck.
Designed with you in mind.

The Full Picture with a 360° View

The FreedomX™ Couch Extensions provide a rigid and lightweight platform for an unobstructed 360° treatment of a patient's head and neck. Our couch extensions feature integrated pitch adjustability which can be combined with the FreedomX™ Roll Adapter (BM02-T) giving simple and precise adjustment for pitch and roll.

The pitch and roll provided by the FreedomX™ Couch Extension combined with yaw from your couch provides a full 6-axis adjustability for patient position correction based on Cone Beam image verification and/or Surface Image Guided Radiation Therapy (SIGRT or SGRT).



Couch Extension Module



The Couch Extension Module is one of the ten highly-designed FreedomX™ Modules.

FDM-T-SRSHN	Freedom™ 6-Axis Linac Couch Extension Bundle for SRS/SRT and Head & Neck Includes: SB02, SB02-S1, BM02-T, SB-HP (couch specific), STA-12-X & optional THS03-MM	1 pkg.
FDM-T-SRS	Freedom™ 6-Axis Linac Couch Extension Bundle for SRS/SRT Includes: SB02-S1, BM02-T, STA-FILL-X, Freedom Spares Kit	1 pkg.
FDM-T-HN	FreedomX™ Treatment Module for Head & Neck Only Includes: SB02, BM02-T, SB-HP (couch specific), STA-12, & optional THS03-MM	1 pkg.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.