




FREEDOMX™

THE FREEDOMX™ TOTAL BODY PRECISION PATIENT IMMOBILIZATION SYSTEM FOR SRS, SBRT AND MORE.

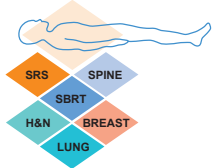

A solid foundation is everything.


The FreedomX™ Comprehensive Overlay is the foundation of the FreedomX™ System and has been engineered with you in mind. At 1.7kg/3.8 lbs, this low density overlay weighs less than a bag of sugar but is as strong as steel.

Compatible with all treatment couches using the corresponding indexing bar, the FreedomX™ Comprehensive Overlay is a true precision centric patient positioning system.



Unique 5-point tension lock mechanism ensures masks are securely fastened to the overlay with easy removal.

Slots allow for attachment of FreedomX Intuition Modules (for SRS and Head & Neck) and FreedomX Shoulder Retractors.

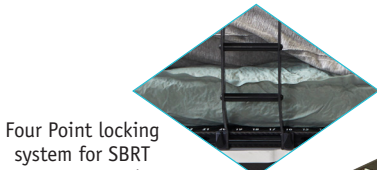

Four Point locking system for SBRT components/modules.




Indexing Bars allow indexing to any couch.



Integrated Table Pad ensures patient comfort. Remove and the X feature will accurately index a posterior VacuumLock™.

Complete compatibility with all robotic couch systems including Elekta™ iGuide Tracking System.

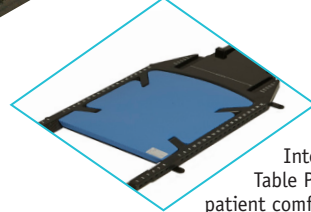

FDM-S-STD-X	FreedomX™ System Comprehensive Overlay Package	1 pkg.
	Includes: LB-SBRT-X, SBRT-FILL-X, CFJ02, couch specific carbon fiber index bars (2)	

LB-SBRT-X	FreedomX™ Comprehensive CT/Treatment Overlay, includes CFJ02	1 ea.
SBRT-FILL-X	FreedomX™ Integrated Table Pads	1 ea.
LB-SBRT-WM-X	FreedomX™ Comprehensive Wall Mount	1 ea.
CFJ02	5-Point Tension Lock Immobilization Clamp	1 ea.
RS02	Shoulder Retractors	1 pair
IB06*	Radio-translucent Indexing Bars	2 ea.


*Note: contact your sales representative to determine the correct part code for indexing bars that match your couch.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™

THE FREEDOMX™ POWERWEDGE MODULE

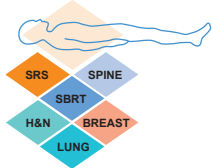

Additional patient comfort with the PowerWEDGE + Knee Positioner
 Designed with you in mind!




Maximizing Patient Comfort

Not all patients can tolerate a flat position. The PowerWEDGE provides enhanced comfort and improves anatomical setup positions with the proven FreedomX™.

As a FreedomX™ module it can be used with virtually any supine arms up or arms down patient positioning and immobilization setup including Breast, Head and Neck, Thorax, Lung, Spine, SRS, SBRT, IGRT, Palliative and COPD patients. The PowerWEDGE is compatible with all imaging modalities and indexes directly to any treatment couch while providing greater comfort for standard treatments to the most complex high precision patient setups. It is light weight and easy to use.

PowerWEDGE Module

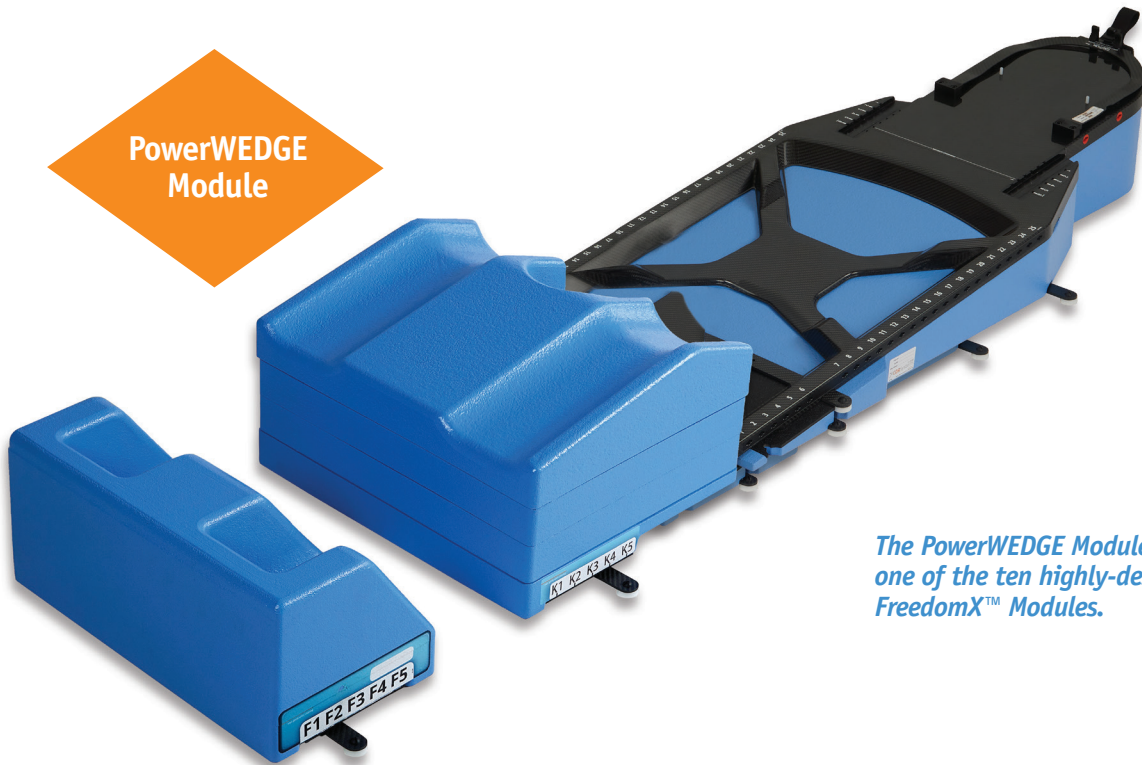

The PowerWEDGE Module is one of the ten highly-designed FreedomX™ Modules.


PWR WDG	FreedomX™ PowerWEDGE	1 pkg.
KN2-WDG-2	FreedomX™ 2 nd Gen Height Adjustable Knee Positioner Indexable without spacer	1 ea.
F2-WDG	FreedomX™ 2 nd Gen Indexable Foot Positioner	1 ea.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™

THE FREEDOMX™ SBRT MODULE

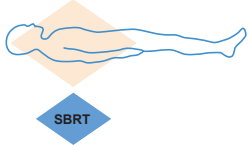


Thoughtful design and engineering allow you to setup with precision and ease.
Designed with you in mind!

Full Body Positioning

As well as Cranial SRS or Head & Neck treatments, the FreedomX™ System can be configured for SBRT type treatment using several additional components that attach and index in conjunction with the FreedomX™ Comprehensive Overlay.

Designed to utilize Vacuum Lock™ cushions, the FreedomX™ System can be manipulated to ensure your patient is extremely comfortable while giving you the highest levels of precision in SBRT patient positioning. The FreedomX™ System is easily configured for Arms Up or Arms Down ensuring your patient will be as comfortable as possible.

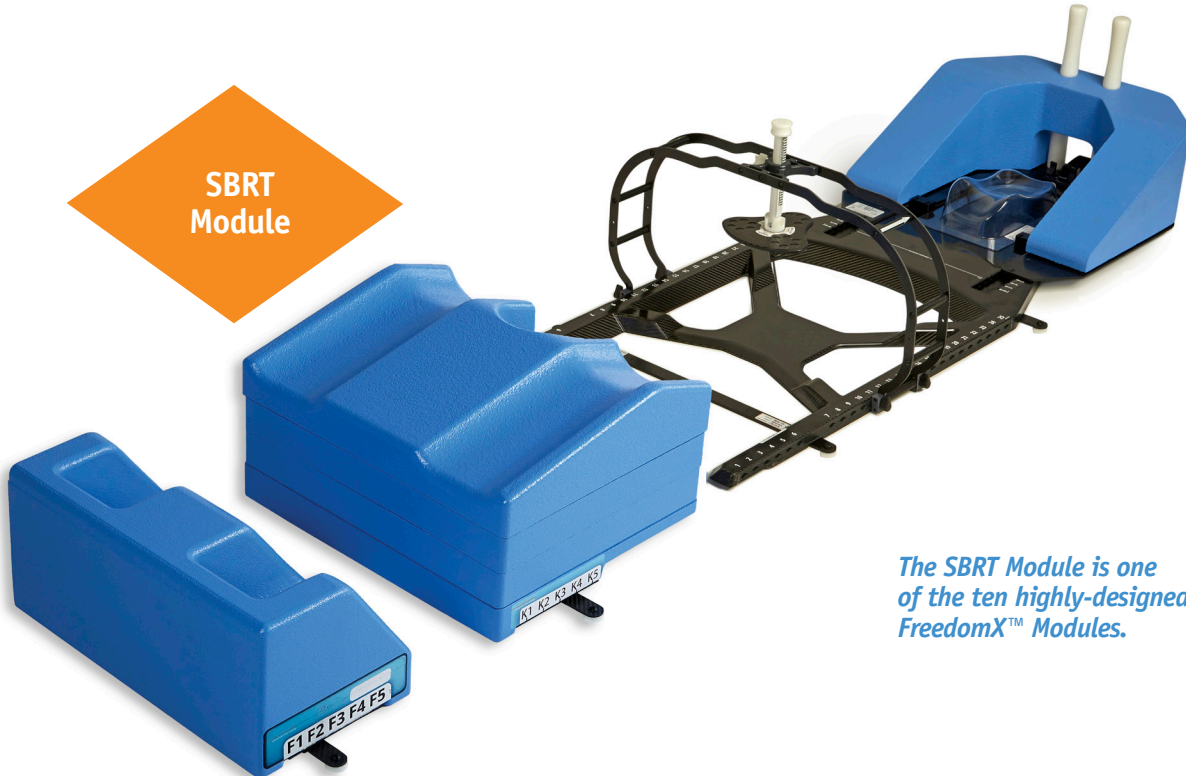

The SBRT Module is one of the ten highly-designed FreedomX™ Modules.

FDM-S-SBRT-X FreedomX™ SBRT Module for Comprehensive Overlay (LB-SBRT-X) 1pkg.
Includes: AU-3, KN2-WDG-2, F2-WDG, SBRT-BS-X, SBRT-BL-X, CMP-S-X, CMP-L-X, IB06 type Indexing Bars x2 (specify couch)


AU-3-SA FreedomX™ Standalone Arms Up module, indexes directly to most couches 1 ea.
Contact your sales representative for couch compatibility.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™


THE FREEDOMX™ SPINE MODULE (REQUIRES SBRT MODULE)

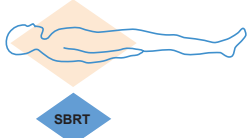

Significantly faster patient setup, planning and treatment.
Designed with you in mind!


Precise Incorporating Straight Alignment

The FreedomX™ Spine-Align Module adds new functionality to our proven low attenuating carbon fiber bridges. With Spine-Align you have the functionality of side compression paddles ensuring your patient is aligned and indexed centrally and reproducibly at simulation and in every treatment. A benefit of using Spine-Align is that any lateral corrections based on an initial scout image will now be much faster to correct. This allows for significantly faster patient setup and unmatched patient immobilization throughout initial patient setup at simulation and for each treatment fraction.

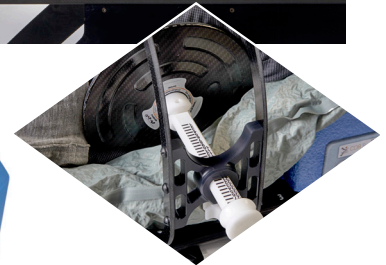

> Spine-Align can be used in any patient treatment position including the arms up position

Spine Module

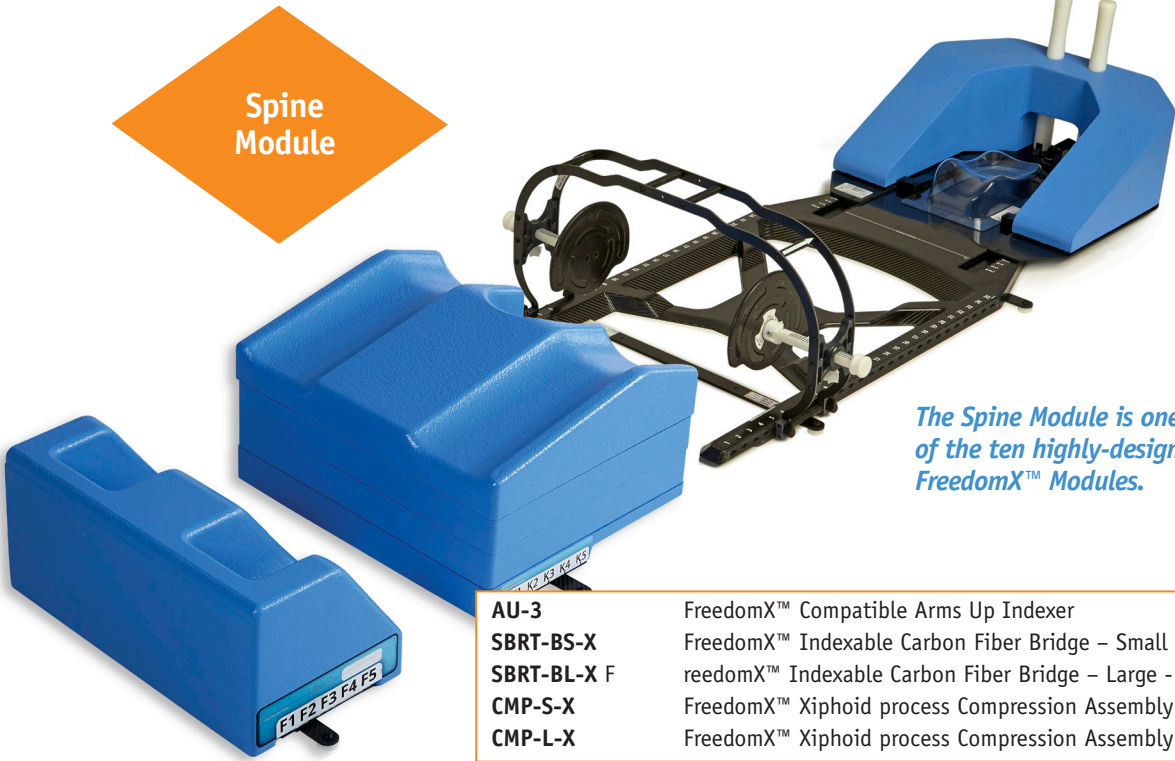

The Spine Module is one of the ten highly-designed FreedomX™ Modules.

AU-3	FreedomX™ Compatible Arms Up Indexer	1 ea.
SBRT-BS-X	FreedomX™ Indexable Carbon Fiber Bridge – Small - fitted for abdominal comp	1 ea.
SBRT-BL-X F	FreedomX™ Indexable Carbon Fiber Bridge – Large - fitted for abdominal comp	1 ea.
CMP-S-X	FreedomX™ Xiphoid process Compression Assembly – Small	1 ea.
CMP-L-X	FreedomX™ Xiphoid process Compression Assembly – Large	1 ea.
SBRT-BS-SP1-X	FreedomX™ Small Indexable Bridge Configured for Patient Spine Alignment	1 ea.
SBRT-BL-SP1-X	FreedomX™ Large Indexable Bridge Configured for Patient Spine Alignment	1 ea.
SP1-M	FreedomX™ Spine-Align Side Positioning Assembly Medium Post	1 ea.
AAS03	Arm/Elbow Adjustable Support Strap	1 ea.
KN2-WDG-2	FreedomX™ 2 nd Gen Height Adjustable Knee Positioner Indexable without spacer	1 ea.
F2-WDG	FreedomX™ 2 nd Gen Indexable Foot Positioner	1 ea.
AU-3-SA	FreedomX™ Standalone Arms Up module, indexes directly to most couches	1 ea.


Contact your sales representative for couch compatibility.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™

THE FREEDOMX™ HAND INDEXER + SHOULDER MODULES

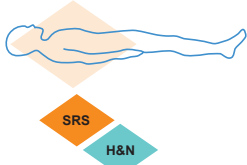

Thoughtful design and engineering allow you to setup with precision and ease. Designed with you in mind!

It's the details that matter

Improve the performance of the most accurate SRS immobilization solutions with the addition of the Hand Indexer and/or the Shoulder-Tech Modules.

Shoulder Tech

> Optimal patient comfort and accuracy of C-spine and shoulders.

Arm/Elbow Support Strap


> Additional stability and comfort is in this one simple solution.

Hand Indexer

> Provides a comfortable position for increased accuracy of arms, shoulder and clavicle.



Hand Indexer + Shoulder-Tech Modules

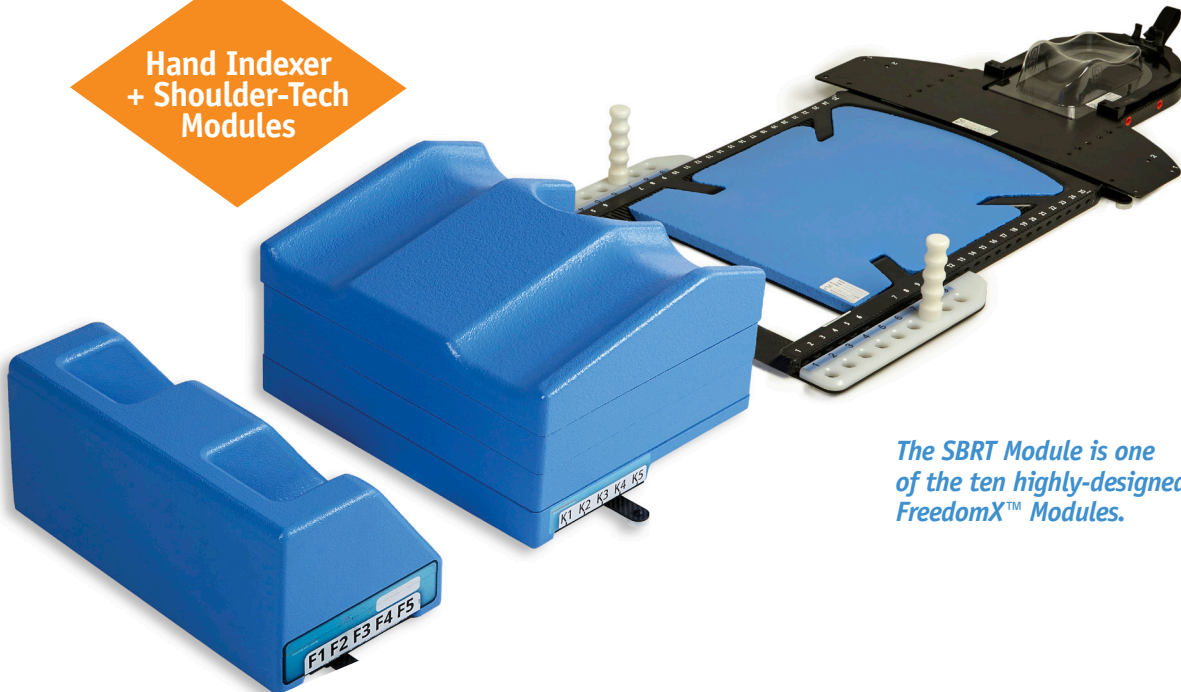

The SBRT Module is one of the ten highly-designed FreedomX™ Modules.


SB-HP	Arms Down Hand Indexer (specify couch)	1 ea.
STA-12-X	Shoulder-Tech – compatible with LB-SBRT-X, LB02, SB02	1 ea.
AAS03	Arm/Elbow Adjustable Support Strap	1 ea.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™

THE FREEDOMX™ DUO™ SRS IMMOBILIZATION SYSTEM

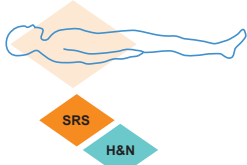

The DUO™ Mask Immobilization System. Ease and accuracy with you in mind.




FitAdjust™ is a quick and easy way to adjust the fit for maximum comfort.

It's the details that matter

CDR System's new FreedomX™ DUO™ SRS immobilization system with FitAdjust™ for instant mask comfort and tightness adjustability, SRS has never been so easy. The FreedomX™ DUO™ SRS System is compatible with most immobilization platforms making upgrading your immobilization as simple as click and go!

The System features:

- > Carbon Fiber Leaf Spring to Ensure Tension
- > Dual Mask Fasteners for Customization
- > Adjustable Latch System to Maximize Comfort
- > Strong and Lightweight Carbon Fiber Base

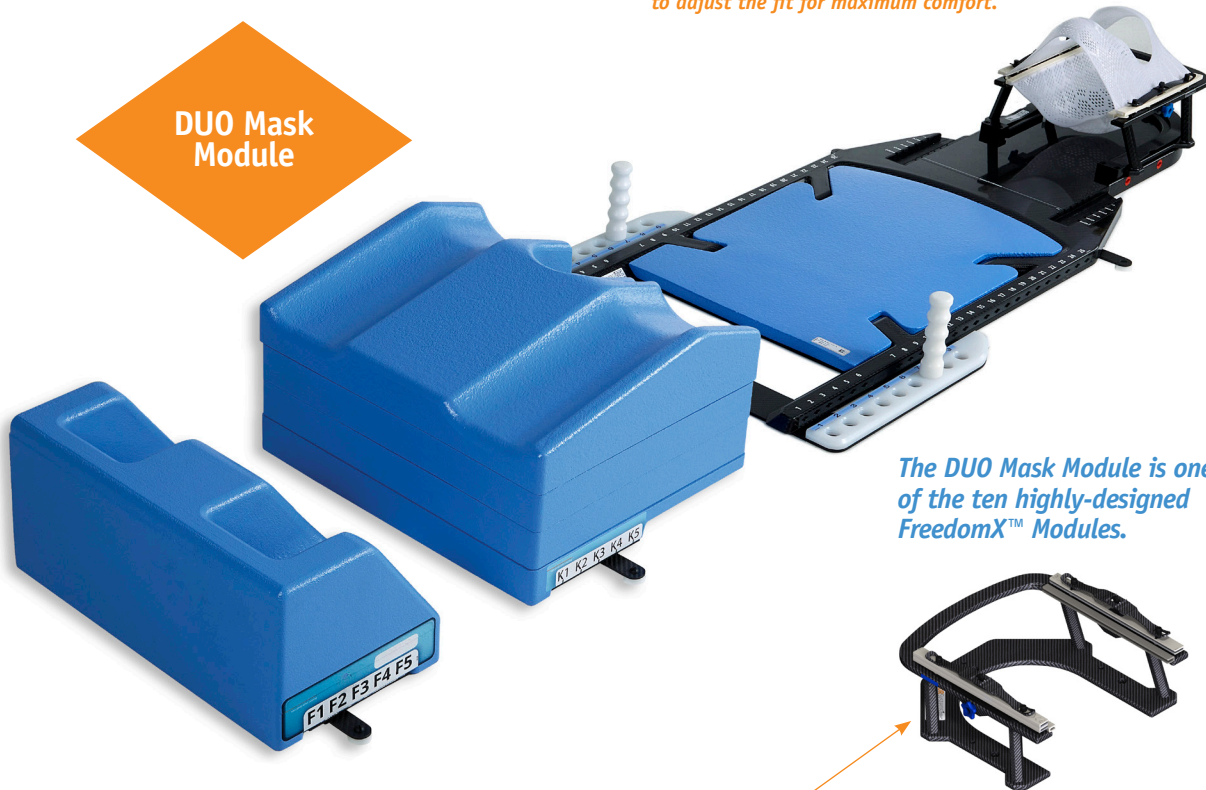

The DUO Mask Module is one of the ten highly-designed FreedomX™ Modules.


DUO SYSTEM:	(Compatible with Type-S™ Overlay)	
DU022-STD	FreedomX™ Duo Mask Head Support Standard	1 pkg.
DU022-LRG	FreedomX™ Duo Mask Head Support Large	1 pkg.
DU022-PED	FreedomX™ Duo Mask Head Support Pediatric	1 pkg.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.




FREEDOMX™


THE FREEDOMX™ DUO SRS IMMOBILIZATION SYSTEM

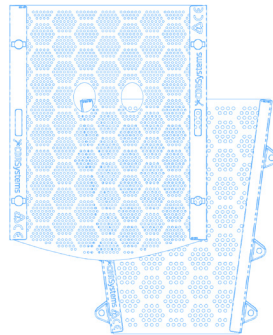

DMS-BASIC

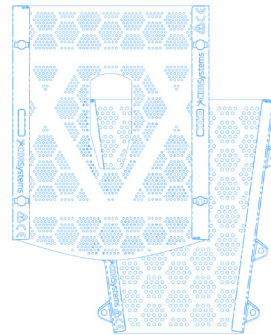

DMS-OPEN

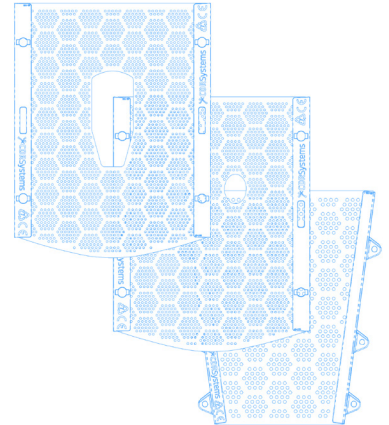

DMS-STER

DMS-BASIC	FreedomX™ Duo Mask Set, Basic Includes 1 Bottom Mask + 1 Basic Top Mask	1 pkg.
DMS-OPEN	FreedomX™ Duo Mask Set, Open Includes 1 Bottom Mask + 1 Open Mask	1 pkg.
DMS-STER	FreedomX™ Duo Mask Set, Stereotactic Includes 1 Bottom Mask + 2 Top Masks	1 pkg.



Shown here: FreedomX™ Duo Mask Set, Open Blue




Shown here: FreedomX™ Duo Mask Set, Open White



FREEDOMX™

THE FREEDOMX™ INTUITION CLICK™ MODULE

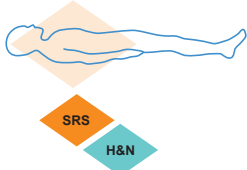

The Intuition Click™ offers superior comfort and accuracy. Ease and accuracy with you in mind.


Precision, Non-invasive and Frameless

The Intuition Click™ provides a system to produce a fully customized, rigid, air equivalent, non-invasive positioning solution that provides full posterior support to a patient's head, C-spine, and shoulders. The completed customized support provides rigidity right up to the coronal plane.

Combined with either a closed mask or open-face mask, the Intuition Click™ has been shown to better the performance of the most accurate SRS immobilization solutions with superior patient comfort.



Intuition Click Module

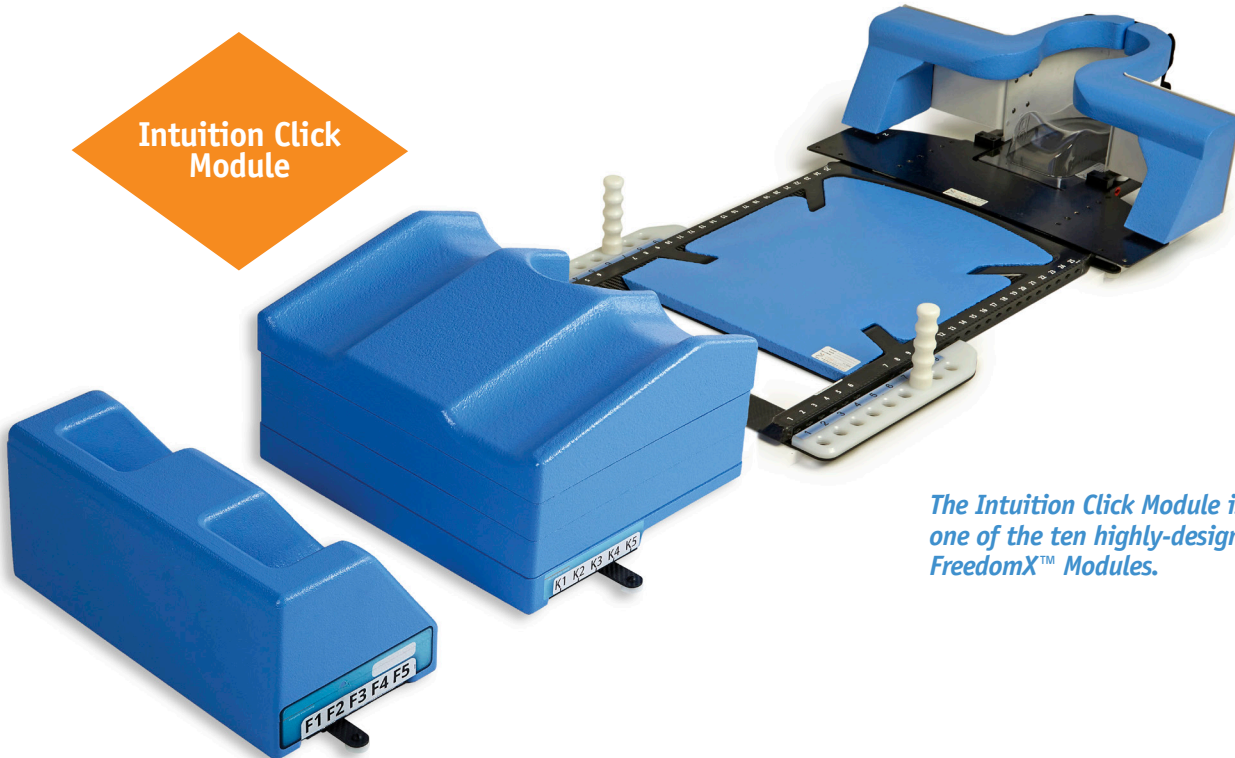

The Intuition Click Module is one of the ten highly-designed FreedomX™ Modules.

The Instaform™ Kit
Easy to mix, the Instaform™ solutions are quickly combined to create one solution for accurate molding.



See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.




FREEDOMX™

THE FREEDOMX™ INTUITION CLICK™ MODULE

The Advantage of an Intuition Click™ created Mold



The Intuition Mold™ creates a true rigid active hold at the cranial stop.

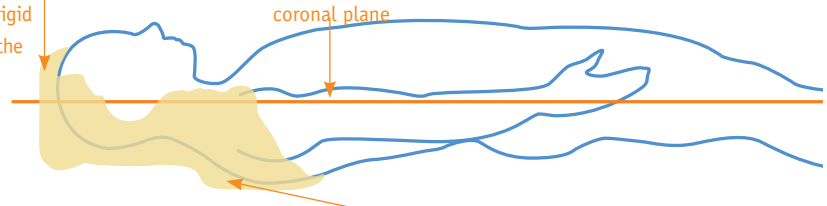

1/ A TRUE MATCH


The Intuition Click™ creates an Intuition Mold™ that is a *true* match to the patient's contour.

2/ AN ACTIVE HOLD AT THE CORONAL PLANE

The Intuition Mold™ is rigid up to the coronal plane for customized support that never relaxes. It creates a true active hold with pressure at the axial part to prevent roll and pitch.

3/ C-SPINE HOLD

The Intuition Mold™ is created using Instaform™. It forms perfectly around the patients head and shoulders – the result is a C-spine hold that tapers to zero.

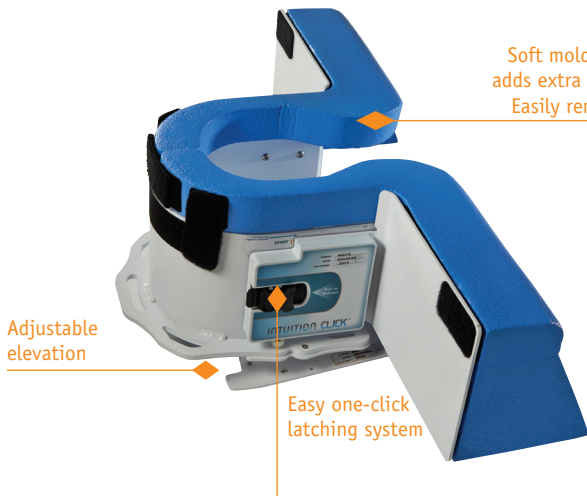


The Ultimate in Accuracy > The Intuition™ Click is a frame used in the patient setup phase to create a custom mold unlike any other!

The Intuition Click Module™ allows you to create a custom Intuition Mold™ that provides full rigid support to a patient's head, neck and shoulders, resulting in superior positioning accuracy.

The creation of the Intuition Click™ Mold using the Instaform™ can be a transformative experience for the patient. Its warmth calms the patient. It has been described as "feeling like a warm bath."

The Intuition Mold™ is as light as air and as solid as the ground we walk on.

FDM-S-SRSHN	Intuition Click™ Module for SRS/SRT + Head & Neck Setup & Simulation	1 pkg.
	Includes: MM14, MM14-IR, LB-HP, STA-FILL-X, STA-12-X, THS03-1-MM	
FDM-S-SRS	Intuition Click™ Module for SRS/SRT & Simulation	1 pkg.
	Includes: MM14-IF, MM14-IR, STA-FILL-X, THS03-1-MM	

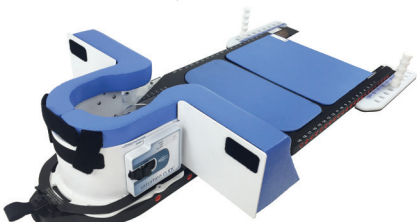


MM14	Intuition Click™ for Custom Head, Neck & Shoulder includes frame + risers	1 ea.
MM14-FI	Intuition Click™ Frame Replacement Foam Insert Kit	1 ea.
MM14-IR	Intuition Click™ Risers – for use with F size tall head supports	1 ea.
STA-12-X	Shoulder-Tech – compatible with LB-SBRT-X, LB02 & SB02	1 ea.
STA-FILL-X	FreedomX™ Cranial SRS Adapters for use with SB02-S1	1 ea.
LB-HP	FreedomX™ Arms Down Hand Indexer – configured to couch	1 ea.
AAS03	Arm/Elbow Adjustable Support Strap	1 ea.



FREEDOMX™

THE FREEDOMX™ MASK STABILIZATION MODULE

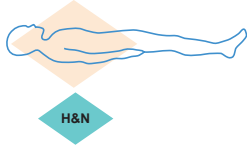

Mask Stabilization offers extra stabilization. Designed with you in mind.


Unparalleled Stabilization

The Mask Stabilization System is an optional add-on that provides additional rigidity and stabilization to meet immobilization protocols. Often used when optimal posterior support is not the focus. This system attaches directly to an immobilization mask and eliminates the need for an invasive bite block system as a means to enhance patient immobilization.

Both the *Wishbone Stabilizer* and *Intuition Stabilizer* index directly to all FreedomX™ Overlays and Couch Extensions using our proprietary LT-Thermoplastics Quick-Attachment process. Both open and closed face masks can be used with FreedomX™ Mask Stabilization System.



Mask Stabilization Module

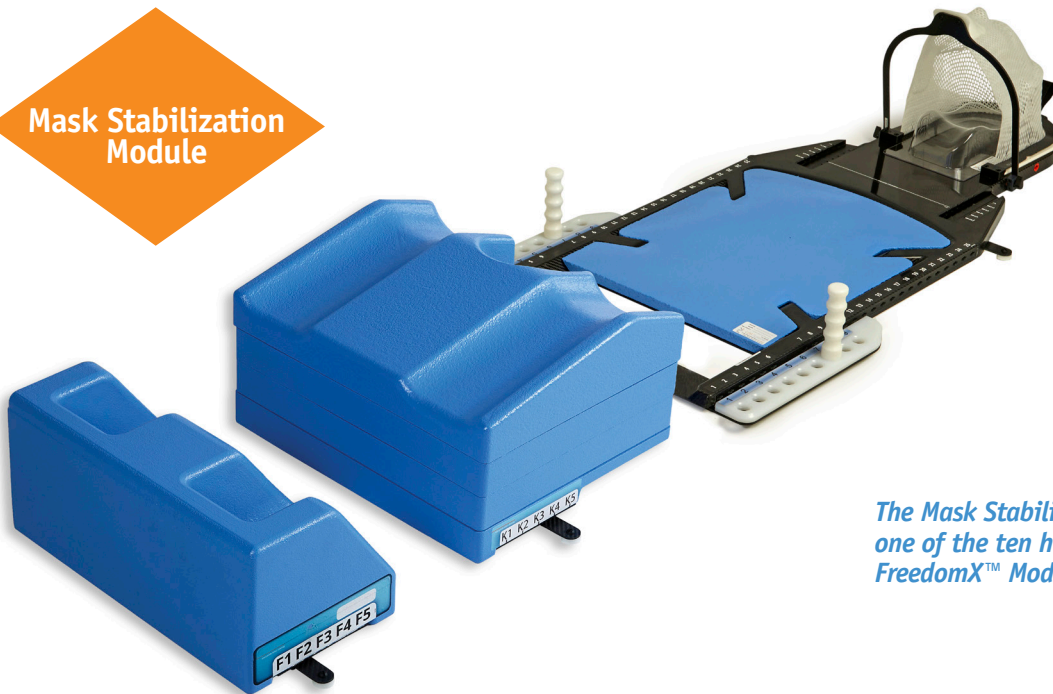

The Mask Stabilization Module is one of the ten highly-designed FreedomX™ Modules.

NWB02-MM	Wishbone Carbon Fiber Stabilizer – large – for use with Intuition	1 ea.
NWB02	Wishbone Mask Stabilizer System – for Freedom Overlays	1 ea.
NWB02-QA	Quick Attachments x 12 ea.	1 pkg.


IFF03	Instaform™ Kit – Standard (4 Kits/Box)	4 ea.
IFF03-B	Instaform™ Kit – Large (4 Kits/Box)	4 ea.
IFF03-D	Instaform™ Kit – Extra Large (8 Kits/Box)	8 ea.
IFF03-E	Instaform™ Kit – Treble Size (6 Kits/Box)	6 ea.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.



FREEDOMX™

THE FREEDOMX™ COUCH EXTENSION MODULE

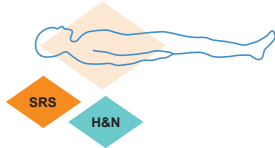

For an unobstructed 360° treatment of a patient's head and neck. Designed with you in mind.

The Full Picture with a 360° View

The FreedomX™ Couch Extensions provide a rigid and lightweight platform for an unobstructed 360° treatment of a patient's head and neck. Our couch extensions feature integrated pitch adjustability which can be combined with the FreedomX™ Roll Adapter (BM02-T) giving simple and precise adjustment for pitch and roll.

The pitch and roll provided by the FreedomX™ Couch Extension combined with yaw from your couch provides a full 6-axis adjustability for patient position correction based on Cone Beam image verification and/or Surface Image Guided Radiation Therapy (SIGRT or SGRT).



Couch Extension Module



The Couch Extension Module is one of the ten highly-designed FreedomX™ Modules.

FDM-T-SRSHN	Freedom™ 6-Axis Linac Couch Extension Bundle for SRS/SRT and Head & Neck Includes: SB02, SB02-S1, BM02-T, SB-HP (couch specific), STA-12-X & optional THS03-MM	1 pkg.
FDM-T-SRS	Freedom™ 6-Axis Linac Couch Extension Bundle for SRS/SRT Includes: SB02-S1, BM02-T, STA-FILL-X, Freedom Spares Kit	1 pkg.
FDM-T-HN	FreedomX™ Treatment Module for Head & Neck Only Includes: SB02, BM02-T, SB-HP (couch specific), STA-12, & optional THS03-MM	1 pkg.

See cdrsys.ca for more information on the FreedomX Modules and FreedomX Accessories.