



S-Type Thermoplastics

For S-Type (Push-Pin) Overlays Immobilization Masks

Our Blue S-Type LT-Thermoplastics™ from CDR Systems supplement our versatile range of masks with both standard perforation and reinforced designs, ensuring secure, precise positioning for all anatomical areas.

Engineered for superior rigidity and reproducibility, these masks support consistent patient positioning throughout treatments. Our Open-Faced Masks are specifically crafted to optimize the accuracy of your Surface Guided Radiation Therapy (SGRT) system, helping to deliver more precise outcomes.



LOW SHRINKAGE



SGRT COMPATIBLE



LOW ACTIVATION
TEMPERATURE



X-RAY TRANSLUCENT



NON-STICK



LATEX FREE



PHTHALATES FREE



MR COMPATIBLE


S-Type LT-Thermoplastic Head Masks

2.4 MM

SH-STYP-24-B	S-Type Head Mask Blue	5 pk
SH-STYP-24R-B	S-Type Head Mask Blue Reinforced	5 pk
SH-STYP-24RO-B	S-Type Head Mask Open Face Blue Reinforced	5 pk


3.2 MM

SH-STYP-32-B	S-Type Head Mask Blue	5 pk
SH-STYP-32R-B	S-Type Head Mask Blue Reinforced	5 pk
SH-STYP-32RO-B	S-Type Head Mask Open Face Blue Reinforced	5 pk

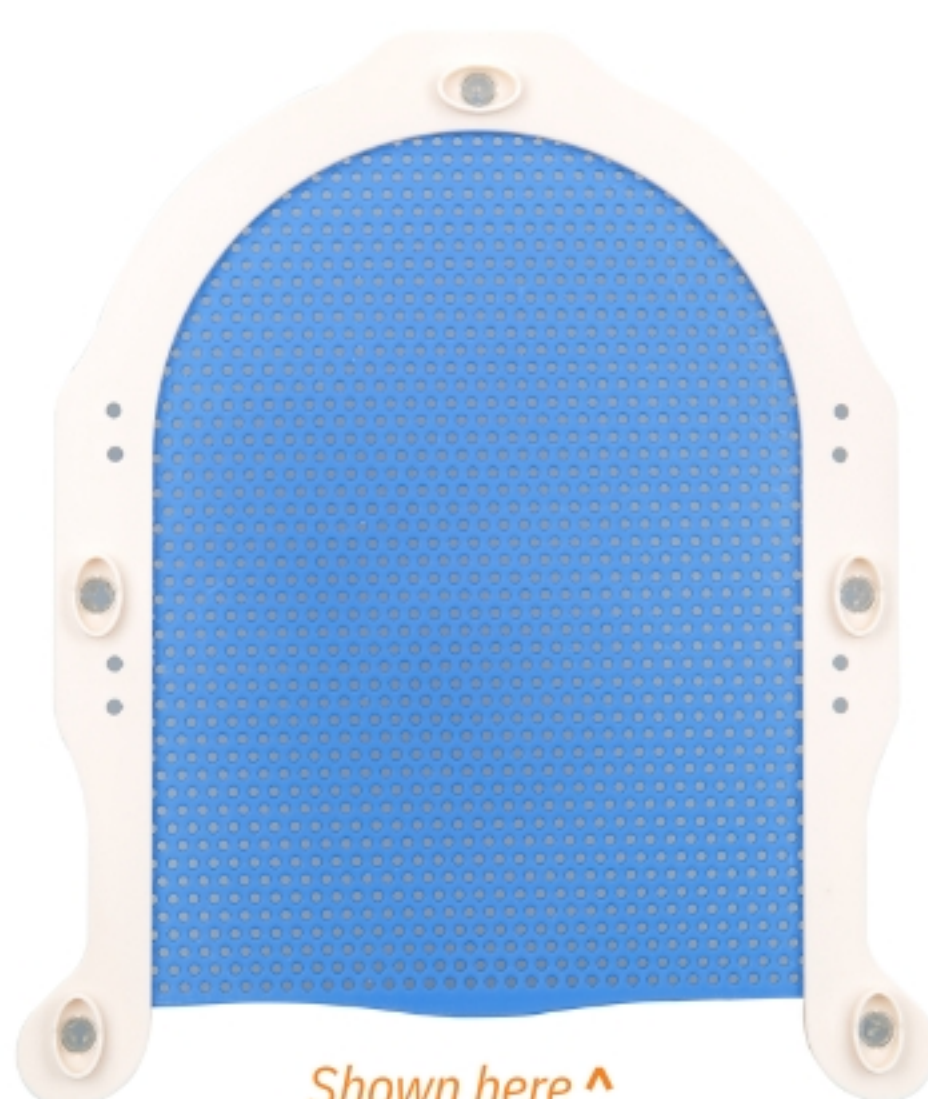


Shown here ^

SH-STYP-24-B
(Available in 2.4 or 3.2mm)

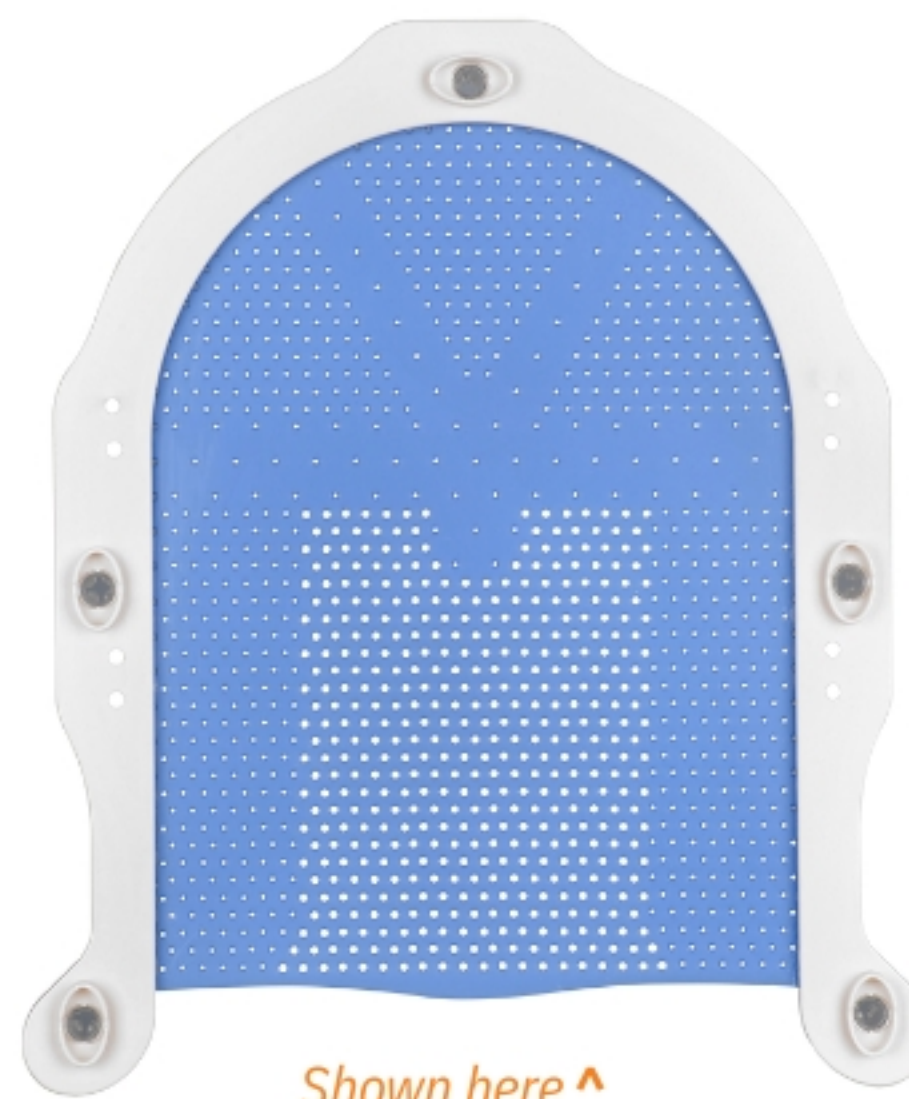

Shown here ^

SH-STYP-24R-B
(Available in 2.4 or 3.2mm)

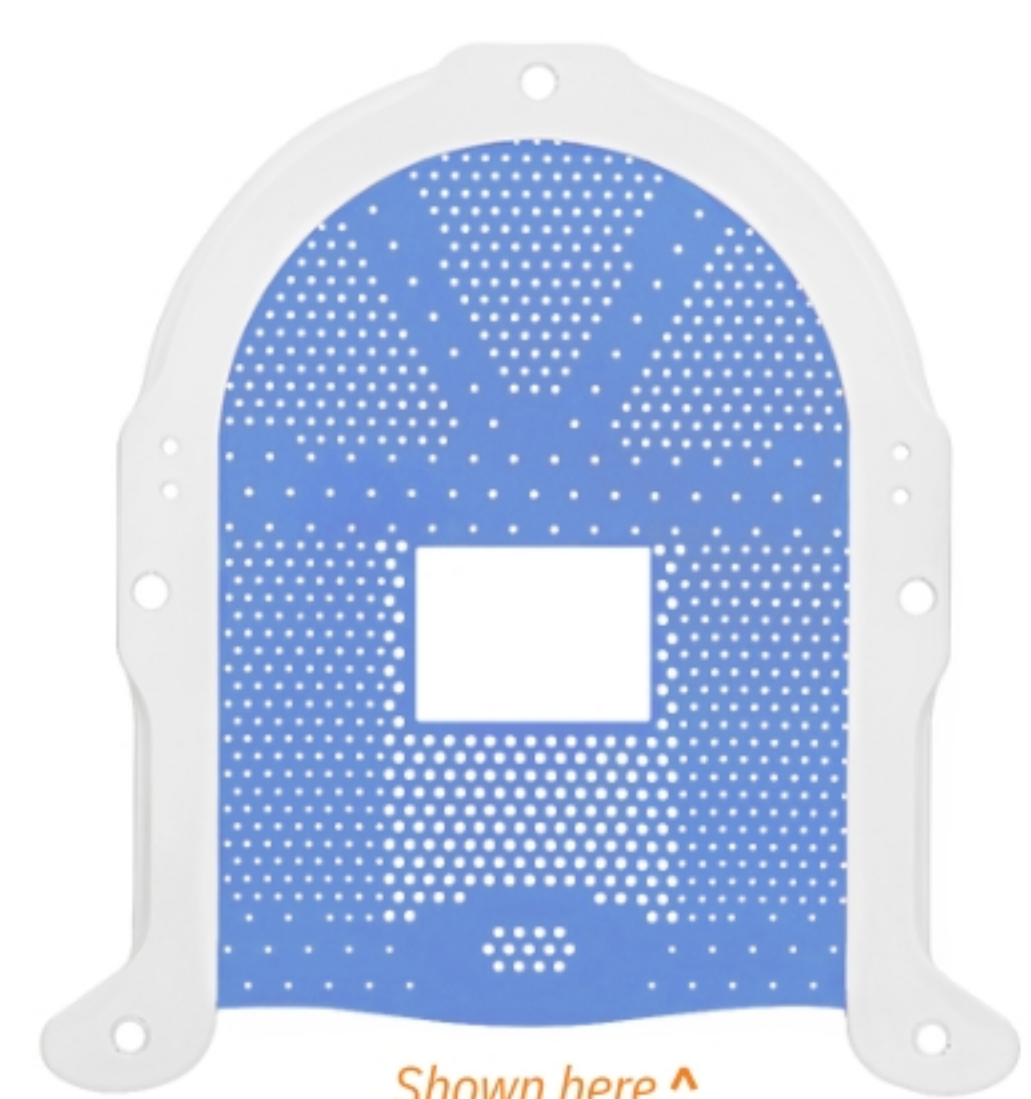

Shown here ^

SH-STYP-24RO-B
(Available in 2.4 or 3.2mm)



S-Type LT-Thermoplastic Head & Shoulder Masks

2.4 MM


ST-STYP-24-B	S-Type H&S Mask Blue 2.4mm	5 pk
ST-STYP-24R-B	S-Type H&S Mask Blue 2.4mm Reinforced	5 pk

2.4 MM Pediatric

ST-STYP-24P-B	S-Type H&S Mask Blue 2.4mm for PEDIATRIC	5 pk
----------------------	-----------------------------------------------------	-------------

3.2 MM


ST-STYP-32-B	S-Type H&S Mask Blue 3.2mm	5 pk
ST-STYP-32R-B	S-Type H&S Mask Blue Reinforced 3.2mm	5 pk
ST-STYP-32RO-B	S-Type H&S Mask Open Face Blue Reinforced 3.2mm	5 pk
ST-STYP-32RON-B	S-Type H&S Mask Open Face Open Neck Blue Reinforced 3.2mm	5 pk

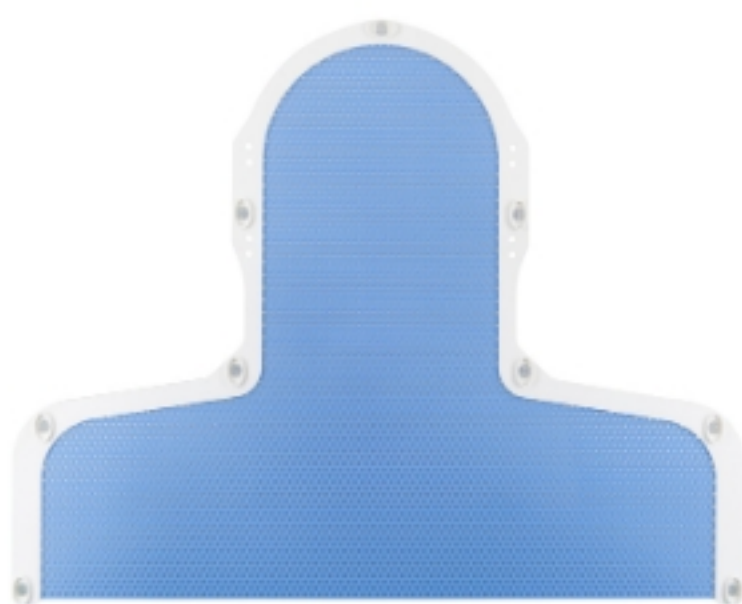


Shown here ^
ST-STYP-24-B
(Available also in 3.2mm)




Shown here ^
ST-STYP-24R-B
(Available also in 3.2mm)

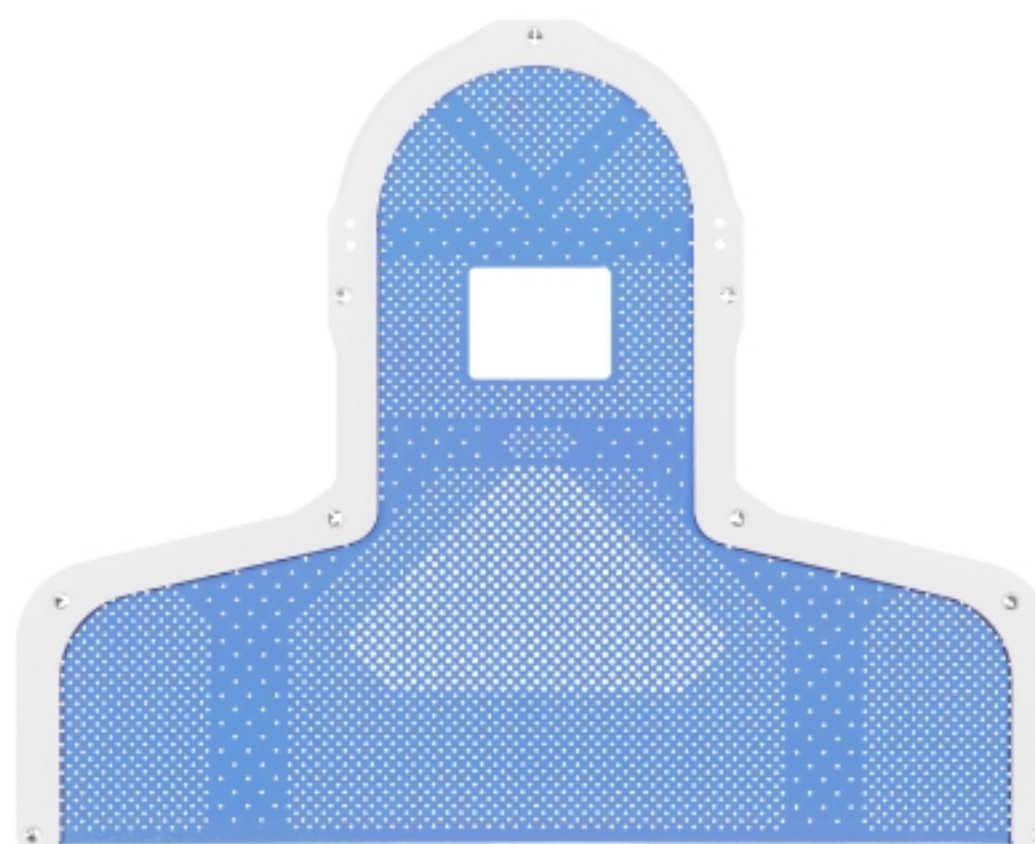

Shown here ^
ST-STYP-32RO-B

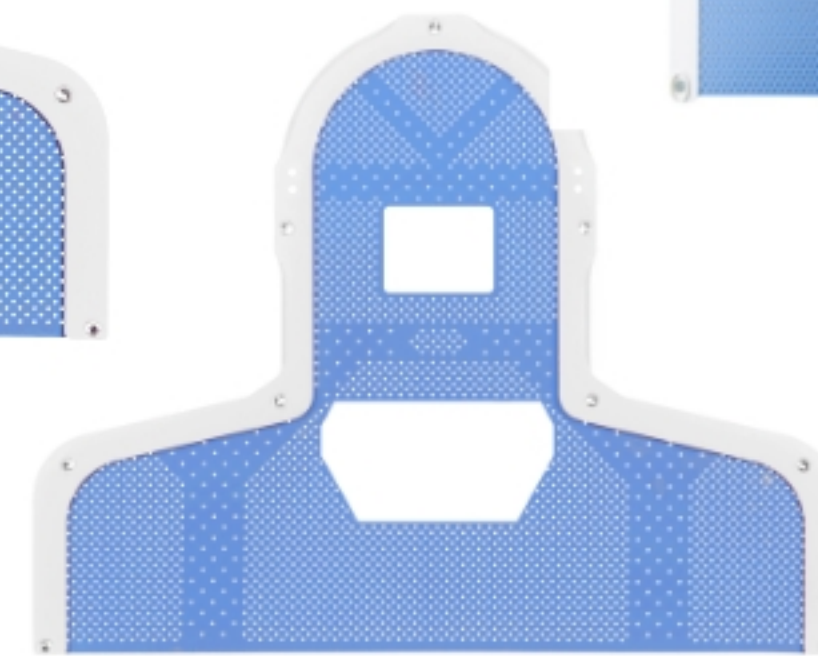

Shown here ^
ST-STYP-32RON-B

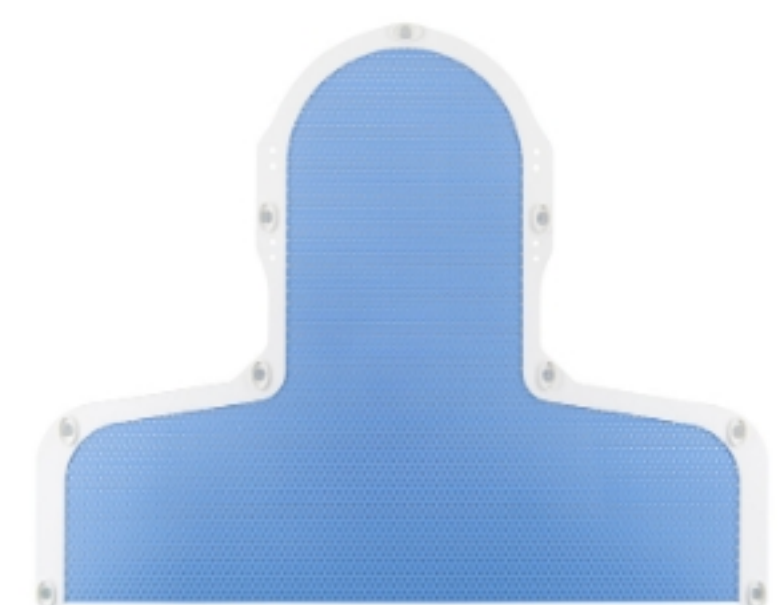

Shown here ^
ST-STYP-24P-B

Moulding & Reference Data

PRODUCT REFERENCE	IDEAL ACTIVATION TEMP	HEATING TIME (Water Bath)	MODELING TIME	HARDENING TIME
LT-Thermoplastic Blue 2.4mm	70C°	2-3 Min	2-2.5 Min	12-15 Min
LT-Thermoplastic Blue 3.2mm	70C°	3-4 Min	2-3 Min	12-15 Min

PRODUCT REFERENCE	DENSITY (g/cm3)	ABSORPTION 6MV	ABSORPTION 10 MV	SHRINKAGE (after 24 hrs)
LT-Thermoplastic Blue 2.4mm	1.6	0.35	0.15	0,02
LT-Thermoplastic Blue 3.2mm	1.6	0.35	0.15	0,02

**Data provided is for reference only. Centers should perform measurements again to verify these results.

Recommended Products



Shown here ^

LBST - S-Type Overlay
Compact Overlay for Type S (Push-Pin) Frames



Shown here ^

DHCO-1 -110 - Digital Dry Heat Oven
(220-240V 50/60Hz also avail.)



Shown here ^

WB-LRG-3 - Digital Waterbath (120V)
(220-240V 50/60Hz also avail.)



Shown here ^

TOTIM - TTM3010M
Custom Patient Posterior Cushion